



OPEN ACCESS

ORIGINAL RESEARCH

Experimental evaluation of stent retrievers' mechanical properties and effectiveness

Paolo Machi,¹ Franck Jourdan,² Dominique Ambard,² Cedric Reynaud,² Kyriakos Lobotesis,³ Mathieu Sanchez,² Alain Bonafé,¹ Vincent Costalat¹¹Centre Hospitalier Universitaire Gui de Chauliac, Montpellier, France²Department of Mechanical Engineering, University of Montpellier, France³Imperial College Healthcare NHS, London, UK**Correspondence to**Dr Paolo Machi, Centre Hospitalier Universitaire Gui de Chauliac, 80 Avenue Augustine Fliche, Montpellier 34000, France; paolo.machi@gmail.com

Received 15 December 2015

Revised 22 February 2016

Accepted 26 February 2016

Published Online First

25 March 2016

ABSTRACT**Background** Five randomized controlled trials recently appeared in the literature demonstrating that early mechanical thrombectomy in patients with acute ischemic stroke is significantly related to an improved outcome. Stent retrievers are accepted as the most effective devices for intracranial thrombectomy.**Objective** To analyze the mechanical properties of stent retrievers, their behavior during retrieval, and interaction with different clots and to identify device features that might correlate with the effectiveness of thrombus removal.**Materials and methods** All stent retrievers available in France up to June 2015 were evaluated by mechanical and functional tests aimed at investigating the variation of their radial force and their behavior during retrieval. Devices were also tested during in vitro thrombectomies using white and red experimental thrombi produced with human blood. Functional tests and in vitro thrombectomies were conducted using a rigid 3D printed vascular model.**Results** Mechanical tests showed a variation in radial force during retrieval for each stent. A constant radial force during retrieval was related to continuous cohesion over the vessel wall and a higher rate of clot removal efficacy. All stent retrievers failed when interacting with white large thrombi (diameter ≥ 6 mm).**Conclusions** None of the tested devices were effective in removing white clots of large diameter (≥ 6 mm). Constant radial force during retrieval allows constant cohesion to the vessel wall and pressure over the clot; such features allow for a higher rate of clot removal.

during retrieval, and thrombus biomechanics and consistencies.

In this study we experimentally evaluated the interaction between STRs and different types of thrombi. Furthermore, we analyzed mechanical properties of various devices in order to investigate their interaction with the clot and their behavior during retrieval. The aim of this study was to try to identify any properties of the mechanical device that correlate with a higher rate of thrombus removal.

MATERIAL AND METHODSSTRs analyzed in this study are listed in [table 1](#). Each device was evaluated by mechanical and functional tests.**Mechanical tests**

Mechanical tests aimed to measure the radial force exerted by the stent in two specific conditions: upon deployment and during the retrieval.

Flat plate compression testThis test measured the STR outward radial force density⁶; it corresponds to the force exerted by the device against a flat surface during compression. Tests were performed at room temperature with a tensile test machine as follows: the STR was placed between two flat plate transducers and the radial force was measured while the two plates compressed the device. For this study, the radial force was measured while the STRs were compressed to 50% of their labeled diameter. Compression tests consisted of five consecutive cycling compressions and were repeated three times for each device. The mean value of the force exerted by the machine to compress the STR at each cycle was used to calculate the radial force (F) per unit length (L) expressed in N/mm as follows: F (average)/L (after compression).**Pull up traction test**The *Pull up traction test* aimed to evaluate how the outward radial pressure exerted by the STRs varied during retrieval along tubes of different diameters. Tests were performed using a tensile test machine as follows: the device was deployed within a silicone tube maintained by a rigid scaffold and the push wire was connected to the traction transducer arm ([figure 1](#)). Silicon tubes were previously filled with saline solution heated to 37°C. The transducer arm retrieved the devices with a velocity of 2 mm/s for a distance of 5 cm. The force F exerted by the**INTRODUCTION**The results of five randomized controlled trials recently published showed that early mechanical thrombectomy in patients presenting with acute ischemic stroke due to a large vessel occlusion is related to improved functional outcome and reduced mortality as compared with standalone intravenous fibrinolysis (IVF).^{1–5} Stent retrievers (STRs) were used mainly in these trials and are nowadays recognized as the most effective devices for intracranial thrombectomy. The STR is delivered via a femoral artery approach through a microcatheter and deployed within the occluding thrombus. Subsequently, the device and the microcatheter are gently recovered as a single unit from the patient's body. The interaction between device and thrombus is influenced by multiple factors: device mechanical characteristics, device behavior

CrossMark

To cite: Machi P, Jourdan F, Ambard D, et al. *J NeuroInterv Surg* 2017;**9**:257–263.

Table 1 List of thrombectomy devices evaluated in this study

Device	Size*
Trevo Provue (Stryker, Kalamazoo, Michigan, USA)†	4–20/3–20
Catch (Balt, Montmorency, France)‡	3–15/4–20/6–30
Eric (Microvention, Aliso Viejo, California, USA)†	3–20/4–24/6–44
Preset (Phenox, Bochum, Germany)‡	4–20/6–30
Preset LT (Phenox, Bochum, Germany)‡	3–20/4–20
Embotrap (Nauravi, Galway, Ireland)†	5–21
Separator 3D (Penumbra Inc, Alameda, California, USA)†	4.5–26
Revive (Codman, Raynham, Massachusetts, USA)†	4.5–22
Mindframe (Medtronic, Irvine, California, USA)†	3–23
Solitair FR (Medtronic, Irvine, California, USA)‡	4–20/6–30


*For each device, the first value refers to the nominal diameter, the second value refers to the length expressed in mm.

†Complete axial section device.

‡Incomplete axial section device.

arm during the retrieval was measured in newtons. Therefore, the average of these values (obtained after 10 tractions) was used to derivate the radial pressure R_N exerted by the stent toward the vessel wall. This test was conducted using tubes of 1.5 and 3.5 mm inner diameter and was repeated 10 times for each STR. The radial pressure exerted by the stent on the silicone pipe was calculated considering the contact pressure. Owing to the friction, this contact pressure was normal to the internal surface of the silicone pipe but also had a tangential component R_T . The total force F measured by the tensile test machine was the sum of tangential contact pressures R_T over the entire internal surface of the silicone pipe. Assuming a uniform contact pressure this can be shown as

$$F = 2\pi rLR_T$$

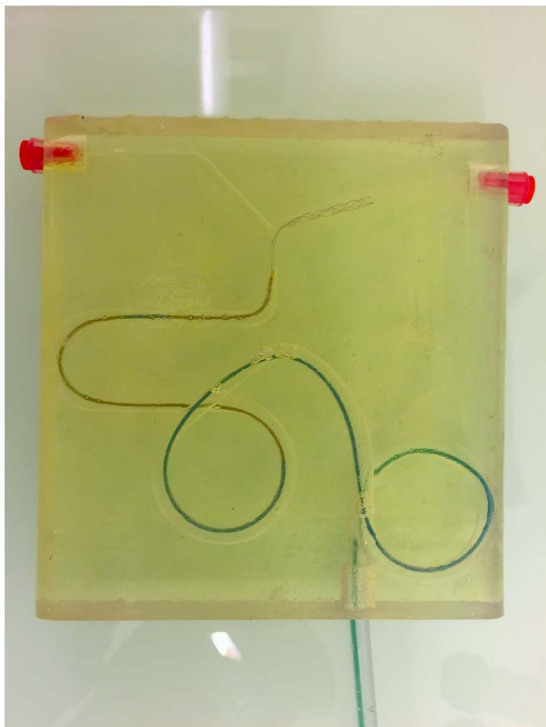**Figure 1** The rigid 3D printed model used for the experiments.

where r is the internal radius of the silicone pipe and L the length of the deployed stent.


According to the Coulomb friction law, there is a relationship between R_N and R_T when the stent is sliding: $R_T = \mu R_N$, which depends on the friction coefficient. Finally, the radial pressure R_N is obtained as follows

$$R_N = F / (2\pi rL\mu)$$

It is expressed in pascals (N/m^2). For our calculations, the friction coefficient was not measured; it was chosen as equal to 0.4, according to values given in the paper by Ranc *et al.*⁷

Functional tests

Functional tests aimed to visually evaluate the ability of an STR to remain in close apposition to the wall of the vessel and to maintain the thrombus engaged within its struts during the retrieval. Tests were performed using a rigid 3D printed vascular model reproducing the brain anterior circulation ([figure 2](#)). Artificial vessels were produced with realistic dimensions but magnified tortuosity in order to evaluate the apposition of the device in unfavorable conditions; the diameter of the middle cerebral artery of the model measured 2.5 mm, the terminal carotid 3 mm, from the carotid syphon to the cervical carotid segment the diameter increased from 4 to 6 mm. During the experiments the model was continuously flushed with saline solution previously heated at 37°C in order to allow the optimal expansion of devices made of nitinol alloy. Two authors with experience in thrombectomy procedures (PM and VC) performed the functional tests, each experiment was filmed and two other authors (FJ and DA) also visually analyzed the results.

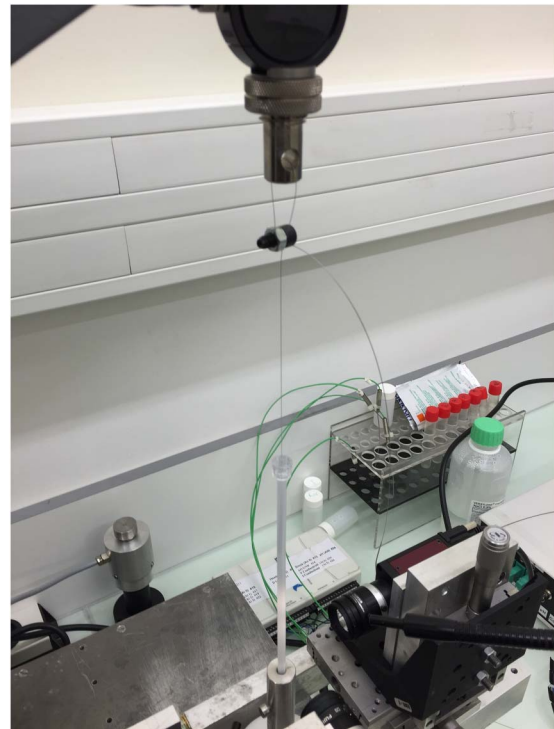**Figure 2** 'Pull up traction test' setting: the figure shows a stent retriever deployed within a silicone tube maintained by a rigid scaffold; the push wire is connected to the traction transducer arm of the tensile test machine.

Retrieving test

The *Retrieving test* was conceived to visually evaluate the ability of an STR' to remain in close apposition to the wall of the vessel during the retrieval. STRs were delivered within the middle cerebral artery of the model. A few minutes after deployment, allowing for their complete expansion, devices were gently retrieved along the vascular model to the guiding catheter connected to the model. The test was repeated 10 times for each stent.

In vitro thrombectomies

This test was aimed to visually evaluate the ability of an STR' to maintain the thrombus engaged within its struts during the retrieval. In order to identify any difference in the behavior of an STR' when interacting with thrombi of different features,⁷ red and white artificial thrombi were employed. In order to ensure that thrombus removal was uniquely due to the device, neither proximal flow arrest nor manual aspiration was performed during retrievals.^{8 9} Artificial thrombi were placed within the middle cerebral artery or terminal carotid artery of the model. Stent retrievers were delivered across the thrombus using a standard technique.¹⁰ The distal extremity of the devices was delivered approximately 1 cm distally to the thrombus. After 3–5 min of embedding time, The STRs were gently retrieved into the guiding catheter. In vitro *thrombectomies* were conducted with red and white artificial thrombi of three diameters: small (2 mm), middle (4 mm), and large (6 mm). Devices were tested five times for each thrombus type; the rate of successful clot removal was calculated singularly for each clot type and as a whole for 'red' or 'white' clots.

Thrombi preparation

Thrombi were generated using human blood obtained from 10 healthy volunteers with no history of use of antiplatelet medication or other medication that might interfere with normal coagulation. Each volunteer provided written informed consent before blood samples were obtained.

White 'stiff' thrombi

Blood samples were collected into tubes containing 3.2% sodium citrate solution (1:9 in volume) used for blood anticoagulation. In order to obtain platelet-rich plasma,¹¹ tubes containing whole blood were centrifuged at 350 g for 10 min at 22°C. After centrifugation, the plasma layer was extracted for establishment of thrombus while the erythrocyte layer was discarded. CaCl₂ was added to the platelet-rich plasma in a ratio of 100 µL:1000 µL to reverse the effect of the sodium citrate. The recalcified plasma was then incubated at 37°C for 1 h. Finally, white thrombi formed from platelet-rich plasma were extracted.

Red 'soft' thrombi

Whole blood samples were collected in standard tubes without sodium citrate solution. After collection, tubes were stored at 37°C for 24 h. Subsequently, red thrombi formed from the whole blood were extracted.¹²

Thrombi flat plate compression test

Flat plate compression tests aimed to evaluate whether stiffness varied between artificial red and white thrombi. Similarly to flat plate tests performed for STRs, thrombi of identical cylindrical shape were placed between two flat plate transducers and the force needed to compress the clots to 50% of their height was

measured. Tests consisted of five non-consecutive compression cycles; the mean value of the force exerted by the machine to compress the clot at each cycle was used to calculate the radial force expressed in mN/mm². White clots were nearly 13 times stiffer than red clots; the force needed to compress white clots was 9 mN/mm², whereas the force needed to compress red clots was 0.7 mN/mm².

Statistics

To investigate whether the variations of the radial pressure exerted by the devices within tubes of diameter 1.5 and 3.5 mm, recorded in *Pull up traction tests*, were statistically significant, a Mann–Whitney test was employed.

The rates of successful clot removal recorded for in vitro *thrombectomies* using different devices for red and white clots were compared with the McNemar test. Differences were considered statistically significant for $p < 0.05$.

RESULTS

Mechanical tests

The results of *Flat plate compression tests* and *Pull up traction tests* are reported in [table 2](#).

During *Flat plate compression* the devices tended to become oval in shape hence discharging radial force toward the axis orthogonal to the plates. This phenomenon was more pronounced for the incomplete section devices (Solitaire FR (Medtronic, Irvine, California, USA), Catch (Balt, Montmorency, France) and Preset (Phenox, Bochum, Germany) in our study) since the free edges of the stent circumferentially overlapped under compression. Consequently, this test was not a reliable comparator for all devices in this study; nevertheless it allowed a separate comparison for complete and incomplete section devices. Among complete section devices, comparing STRs of similar diameter, Eric 4–24 (Microvention, Aliso Viejo, California, USA) had the highest radial force followed by Trevo PV 4–20 (Stryker, Kalamazoo, Michigan, USA), Revive 4.5–22 (Codman, Raynham, Massachusetts, USA), Separator 3D 4.5–26 (Penumbra Inc, Alameda, California, USA), and Embotrap 5–21 (Nauravi, Galway, Ireland) (these last two device showed similar radial force). Among incomplete section devices (Solitaire FR 4–20, Catch 4–20, and Preset 4–20) the value of the radial force density was comparable and it was inferior to the lowest value recorded for complete section devices.

The *Pull up traction tests* permitted the evaluation of the radial pressure variation when the STRs are retrieved along tubes of increasing diameter. Among the complete section devices, Eric 4–24/3–20, Trevo PV 4–20/3–20, Mindframe 3–23 (Medtronic, Irvine, California, USA) and Separator 3D 4.5–26 showed significant reduction of the radial pressure when retrieved within tubes of 3.5 mm in comparison with the pressure recorded in tubes of 1.5 mm inner diameter. Similar behavior was recorded for the Catch 3–15/4–20 and Preset LT 3–20, among incomplete section devices ([table 2](#)). Only Preset LT 4–20 and Preset 6–30 showed no significant difference in radial pressure shift when retrieved in tubes of 1.5 and 3.5 inner diameter; $p = 0.77$ and $p = 0.25$, respectively.

Functional tests

The *Retrieving tests* allowed evaluation of the ability of an STR' to remain in close apposition to the wall of the vessel during the retrieval ([table 3](#)). For the purpose of this test, devices were subdivided into (1) devices that remained in close apposition to the vessel wall during the entire retrieval; (2) devices that remained in apposition to the vessel wall but showed a degree of

Ischemic stroke

Table 2 Result of *Flat plate tests* and *Pull up traction tests*

Size	<i>Flat plate tests</i> Radial force density N/mm	<i>Pull up traction tests</i> Radial pressure Pa (N/m ²) (1.5 mm tube)	<i>Pull up traction tests</i> Radial pressure Pa (N/m ²) (3.5 mm tube)
Trevo PV			
4–20	0.01480	920	50* (p=0.050)
3–20	0.00600	300	60* (p=0.030)
Eric			
3–20	0.01100	220	–
4–24	0.01850	1540	340* (p=0.008)
6–44	0.01304	770	400* (p=0.010)
Embotrap			
5–21	0.00642	1430	770* (p=0.014)
Separator 3D			
4.5–26	0.00791	1360	400* (p=0.006)
Revive			
4.5–22	0.01269	1360	850* (p=0.021)
Mindframe			
3–23	0.00451	1250	330* (p=0.034)
Solitaire FR			
4–20	0.00448	1110	530* (p=0.014)
6–30	0.00351	1060	580* (p=0.020)
Preset			
4–20	0.00521	1090	730* (p=0.018)
6–30	0.00368	630	600* (p=0.25)
Preset LT			
3–20	0.00320	1060	160* (p=0.03)
4–20	0.00370	460	480* (p=0.77)
Catch			
3–15	0.00350	840	220* (p=0.049)
4–20	0.00444	1810	360* (p=0.010)
6–30	0.00368	1870	900* (p=0.030)

*p Value for the comparison between force exerted in tubes of 1.5 and 3.5 mm inner diameter during *Pull up traction tests*. Note results obtained with Preset 6–30 and Preset LT 4–20 for which there is not significant shift of the radial pressure when retrieved in tubes of different diameter; p=0.25 and p=0.77, respectively.

elongation around acute vessel angles ([figure 3](#)); (3) devices that completely lost contact with the vessel wall during the retrieval. The following devices remained in close apposition to the vessel

wall during the entire retrieval: Preset 6–30, Catch 6–30, and Solitaire FR 4–20/6–30. The following devices remained in apposition to the vessel wall but demonstrated elongation when retrieved across acute angles: Embotrap 5–21, Separator 3D 4.5–26, Revive 4.5–22, and Eric 4–24/6–44.

The following devices completely lost apposition to the vessel wall during the retrieval: Preset LT 3–20, Eric 3–20, and Mindframe 3–23.

In vitro *thrombectomies* allowed an appreciation of the interaction between devices and thrombi of different consistency (red or white) and sizes (small, moderate, or large), ([table 4](#)).

With large white thrombi (6 mm), the delivery catheter navigated across the clot by passing between the clot and the vessel wall. Once delivered, stent retrievers remained constrained between the clot and the vessel wall ([figure 4](#)). Thereafter, the devices did not expand but slid over the thrombus; only in a few cases was the device effective in a minimal clot displacement. Interacting with a small or moderate size white thrombus (2–4 mm), stents partially expanded after deployment. Nevertheless, as noted in the larger white thrombi, the clot was not penetrated by the stent struts but remained on the side of the devices compressed against the vessel wall. Contrary to larger thrombi, small and middle thrombi were partially engaged by the stent struts and displaced during the retrieval. Nevertheless, the clots did not remain engaged during the entire retrieval but intermittently rolled between the stent and the vessel wall toward the tip of the stent. The effectiveness of the

Table 3 Results of *Retrieving tests*

Device	Size	Close apposition	Elongation	Loss of contact
Trevo PV	4–20		+	
	3–20			+
Eric	3–20			+
	4–24		+	
	6–44		+	
Embotrap	5–21		+	
Separator 3D	4.5–26		+	
Revive	4.5–22		+	
Mindframe	3–23			+
Solitaire FR	4–20	+		
	6–30	+		
Preset	4–20	+		
	6–30	+		
Preset LT	3–20			+
	4–20		+	
Catch	3–15			+
	4–20		+	
	6–30	+		

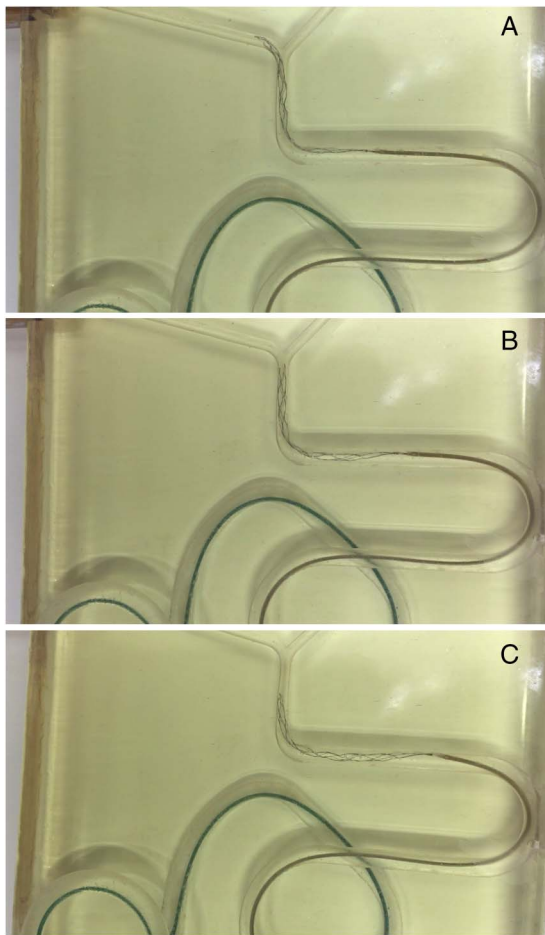


Figure 3 Some devices showed a greater degree of elongation around acute vessel angles: interaction of the proximal, middle, and distal thirds of Solitaire 4–20 (A–C).

removal was secondary to the clot being maintained by the stent struts and the clot rolled between the stent and the vessel wall. The results of *in vitro thrombectomies* performed using white thrombi showed a significant difference between the performance of Preset 6–30 and those of Catch 3–15, Preset LT 3–20, Eric 3–20, and Trevo PV 3–20 (60% vs 20%; $p=0.031$).

A completely different interaction with STRs was noted when *in vitro thrombectomies* were performed with red thrombi. In contrast to white thrombi, the delivery microcatheter navigated through the clot and not alongside it. Furthermore, once deployed, stent retrievers easily penetrated the thrombus with their struts and very often the clot entered within their lumen. Red clots mostly remained engaged with the devices during retrieval, with a higher rate of complete clot removal than recorded for white thrombi. Nevertheless, we noted some differences in embolic complications, which occurred because red thrombi had a tendency to fragmentize during retrievals. Preset 6–30 was effective in red clot removal in 100% of the tests. This rate was significantly higher than the rates obtained with Separator 3D (100% vs 46%; $p=0.008$), Eric 3–20, Preset LT 3–20, Catch 3–15 (100% vs 53%; $p=0.016$), Trevo PV 3–20 and Eric 4–24 (100% vs 60%; $p=0.031$).

Solitaire FR 4–20 and 6–30, Preset 4–20, and Catch 6–30 were effective in red clot removal in 93% of the tests. This rate was significantly higher than the rates of Separator 3D, Eric 3–20, Preset LT 3–20, and Catch 3–15 (93% vs 46%; $p=0.016$ or 93% vs 53%; $p=0.031$).

DISCUSSION

Our study allowed the evaluation of mechanical characteristics of different STRs and permitted a better understanding of how they interact with different types of thrombi.

To date, only a few experimental studies investigating the interaction of thrombectomy devices with artificial thrombi have been carried out.^{13–18} In addition to the device–clot interaction, these studies conducted a comparative evaluation of clot removal efficacy of different thrombectomy devices. In contrast to our study, none of the previous studies evaluated the variation in STR radial pressure during retrieval. Furthermore, these studies predominantly evaluated the interaction of devices with one type of thrombus.

Wenger *et al*¹⁵ evaluated the interaction of four devices (Solitaire, Aperio (Acandis, Pforzheim, Germany), Merci X6, and Merci L5 (Concentric medical, Montain View, California, USA) with artificial clots. The authors found that the migration of the stent struts into the clot is promoted by a strong radial force of the device and by large gaps between filaments. Nevertheless, in our study we found that when delivered within small vessels, across large white clots, devices remain constrained and gaps between filaments did not open. Moreover, devices presenting large gaps between filaments, especially open-cell devices, were more prone to lose contact with the vessel wall when retrieved along sharp angles, disengaging the clot. The authors found that the use of longer STRs improved stable clot engagement. Our results concord with this finding, with the larger devices showing a higher rate of complete clot removal.

Madjidyar *et al*¹⁶ experimentally evaluated device–clot interaction and the influence of distal aspiration on performances of four thrombectomy devices: Solitaire FR, Trevo, Separator 3D, and Aperio. The authors performed experimental thrombectomies, with or without additional distal aspiration, using an artificial vascular model reproducing the human brain anterior circulation. Distal aspiration was performed during retrievals via an intermediate catheter placed at the origin of the middle cerebral artery. The authors found that after stent deployment, the artificial thrombus was pushed against the opposite vessel wall and no integration of the clot into the inner stent lumen could be observed in any device. Moreover, when the stent was pulled back, the clot rolled between the device and the vessel wall and wandered toward the tip of the device. These findings correspond to our observations when analyzing the interaction between devices and small and moderate-sized white clots. Interestingly, the authors noted that when distal aspiration was added to the retrieval process, clot removal efficacy was comparable for all devices. Such results correlate with our observations in this study.

In *pull up traction tests* conducted in 1.5 mm inner diameter tubes we demonstrated comparable radial pressure (>900 Pa) for the majority of the 4 mm diameter devices. Hence, it can be reasoned that all devices have equal clot engagement in the early phase of the retrieval. This might explain why no differences in clot removal efficacy were found when distal aspiration was associated with experimental thrombectomies. On the contrary, we observed that when thrombectomies are performed without any additional aspiration, each device has a different performance, which is related to its own mechanical characteristics. Mechanical tests gave an overview of each STR radial pressure. *Pull up traction tests* performed in tubes of 1.5 mm inner diameter found similar radial pressure for devices of similar diameter (4–4.5 mm). Results of tests performed in tubes of 3.5 mm inner diameter showed different behavior for each

Ischemic stroke

Table 4 Results of in vitro thrombectomies

Device	Size	Large white clot	Medium white clot	Small white clot	%	Large red clot	Medium red clot	Small red clot	%
Trevor PV	4–20	0/5	2/5	3/5	33	3/5	4/5	5/5	80
	3–20	0/5	1/5	2/5	20	3/5	3/5	3/5	60
Eric	3–20	0/5	1/5	2/5	20	3/5	3/5	2/5	53
	4–24	0/5	2/5	3/5	33	2/5	3/5	4/5	60
Embotrap	6–44	0/5	3/5	3/5	40	2/5	5/5	4/5	73
	5–21	1/5*	2/5	3/5	40	3/5	5/5	4/5	80
Separator 3D	4.5–26	0/5	2/5	2/5	26	2/5	2/5	3/5	46
Revive	4.5–22	1/5*	2/5	3/5	40	3/5	5/5	4/5	80
Mindframe	3–23	0/5	2/5	2/5	26	3/5	4/5	3/5	66
Solitaire FR	4–20	0/5	3/5	3/5	40	4/5	5/5	5/5	93‡
	6–30	1/5*	3/5	4/5	53	4/5	5/5	5/5	93‡
Preset	4–20	0/5	3/5	4/5	46	4/5	5/5	5/5	93‡
	6–30	2/5*	3/5	4/5	60†	5/5	5/5	5/5	100§
Preset LT	3–20	0/5	1/5	2/5	20	3/5	2/5	3/5	53
	4–20	0/5	2/5	2/5	26	3/5	3/5	3/5	60
Catch	3–15	0/5	1/5	2/5	20	2/5	3/5	3/5	53
	4–20	1/5*	1/5	2/5	26	3/5	3/5	3/5	60
	6–30	1/5*	3/5	4/5	53	5/5	5/5	4/5	93‡

*Minimal clot displacement.

†The rate of white clot removal found for Preset 6–30 was significantly higher than the rates for Catch 3–15, Preset LT 3–20, Eric 3–20, and Trevor 3–20 (60% vs 20%; $p=0.031$).

‡The rate of red clot removal found for Solitaire 4–20/6–30, Preset 4–20, and Catch 6–30 was significantly higher than the rates for Separator 3D, Eric 3–20, Preset LT 3–20, and Catch 3–15 (93% vs 46%; $p=0.016$ or 93% vs 53%; $p=0.031$).

§The rate of red clot removal found for Preset 6–30 was significantly higher than the rates for Separator 3D (100% vs 46%; $p=0.008$), Eric 3–20, Preset LT 3–20, Catch 3–15 (100% vs 53%; $p=0.016$), Trevor 3–20, and Eric 4–24 (100% vs 60%; $p=0.031$).

device. In this case, most of the tested devices showed a definite decrease of the radial pressure in comparison with values recorded in tubes of 1.5 mm. On the contrary, some incomplete section devices showed constant radial pressure in both 1.5 mm and 3.5 mm inner diameter tubes (Preset 4–20 LT and 6–30, $p=0.77$ and $p=0.25$, respectively, in table 2). Interestingly, the *Retrieving tests* did not show any loss of apposition to the vessel wall during retrieval for such devices. Moreover, a relatively high rate of contact loss was recorded for those devices for which *Pull up traction tests* have shown a sharp decrease of the radial pressure when tested in 3.5 mm inner diameter tubes.

Analyzing the results of in vitro thrombectomies, we noted that devices did not expand when interacting with large white clots (6 mm) but remained constrained between the clot and the vessel wall, and thus during the retrieval, devices slide over the clot without capturing it. None of the devices tested in our study could penetrate and remove large white thrombi (6 mm). Nevertheless, devices could expand and penetrate smaller white thrombi (2–4 mm). During the retrieval, the contact with these clots was not homogeneous but there was a succession of phases in which the clot was engaged by the stent struts and phases in which the clot was disengaged and rolled between the device and the vessel wall.

A completely different interaction was appreciated when devices were tested with red thrombi; in this instance all thrombectomy devices could easily penetrate and incorporate the clot inside their lumen. During the retrieval red clots showed a tendency to fragment and this was related to the occurrence of a higher rate of embolization of distal fragments.

Investigating the interaction devices and thrombi, we found that removal efficacy was related to the device's ability to maintain a constant radial pressure and hence allow gradual expansion and constant apposition during retrieval (eg, Preset 6–30). This reduced the occurrence of distal fragment embolization during the interaction with both white and red thrombi. Such

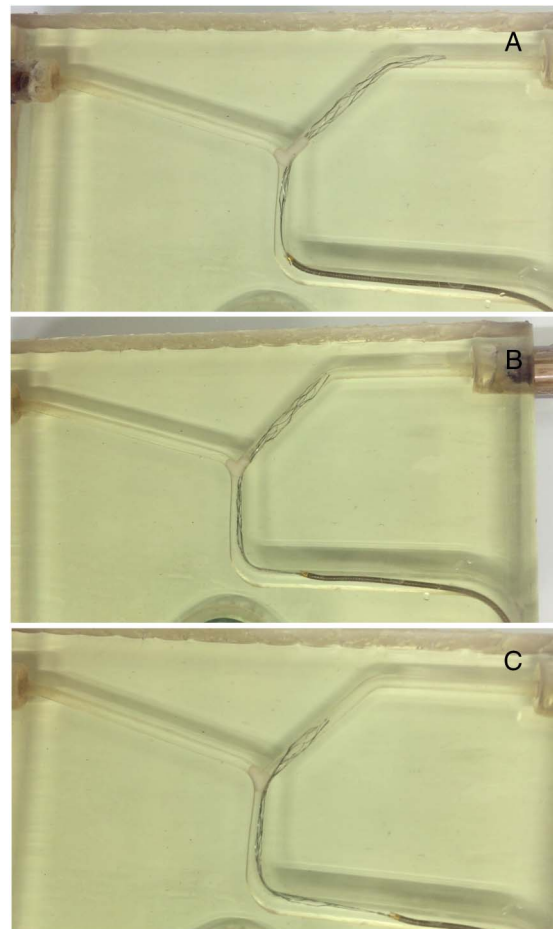


Figure 4 Interaction between Solitaire 4–20 and a large white clot: interaction of the proximal, middle, and distal thirds of the device (A–C).

effective behavior was not recorded for devices for which the radial force reduced sharply during retrieval.

Retrieving tests showed that overall devices flattened when facing sharp vessel angles (90° in our model) (figure 3). Devices that after interaction with sharp angles presented a prompt and complete reopening of the proximal and middle thirds showed better cohesion to the vessel wall and stability during retrieval. Such favorable behavior was mainly recorded for incomplete section devices and, in particular, for devices of larger diameters: Preset 6–30, Solitaire FR 6–30, and Catch 6–30. On the other hand, devices that did not promptly re-expand after interaction with sharp angles lost contact with the vessel wall and jumped abruptly toward the downstream vessels. In our experiments this behavior, mainly recorded for complete section devices, was related to whole clot disengagement and fragment embolization.

Our study has several limitations. First, we used a rigid plastic model and hence during retrievals across the curves and tortuosity of a model, the vasculature did not straighten. Nevertheless, we planned to employ a rigid model in order to exacerbate device performance and to assure identical conditions for all devices. Moreover, the thrombus interaction with the inner surface of the plastic vessel is not comparable with the interaction with the human endothelium, consequently the resistance recorded in our experiments during retrieval was lower than in real vessels. In addition, despite of our best efforts to produce clots with similar features, thrombi slightly varied in size or texture.

CONCLUSIONS

Our results add some innovative information to the current understanding of mechanical thrombectomy.

1. All STRs slide over large white thrombi, with the clot failing to be retrieved. To our knowledge, this thrombus removal failure mechanism has not been reported so far.
2. White, small and medium thrombi are not permanently engaged by the stent struts but roll between it and the vessel wall during retrieval. On the other hand, red thrombi are easily caught by devices but have a tendency to fragment.
3. Constant radial force during retrieval allows constant apposition to the vessel wall and pressure over the clot. Such features allow for a higher rate of clot removal efficacy.
4. Larger devices (6 mm in diameter) have better clot removal performance and increase the probability of capturing the clot.
5. The results of our study have to be compared with further experimental and clinical investigations.

Acknowledgements The authors gratefully acknowledge Professor Thierry Metens for assistance with statistical analysis.

Contributors All authors contributed to planning the research and conducting the experiments. PM and FJ edited the manuscript.

Funding This research received no specific grant from any funding agency in the public, commercial, or not-for-profit sectors.

Competing interests AB unrelated: consultancy: Covidien (consultant for ev3), Stryker; grants/grants pending: Covidien. VC: consultancy: Balt, Codman Neuro-DePuy Synthes, Stryker, MicroVent; payment for lectures (including service on speakers bureaus): Stryker, Balt; payment for development of educational presentations: Covidien.

Provenance and peer review Not commissioned; externally peer reviewed.

Data sharing statement We are agreeable to sharing unpublished data with readers if requested.

Open Access This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/4.0/>

REFERENCES

- 1 Berkhemer OA, Fransen PSS, Beumer D, *et al*. A randomized trial of intraarterial treatment for acute ischemic stroke. *N Engl J Med* 2015;372:11–20.
- 2 Goyal M, Demchuk AM, Menon BK, *et al*. Randomized assessment of rapid endovascular treatment of ischemic stroke. *N Engl J Med* 2015;372:1019–30.
- 3 Campbell BC, Mitchell PJ, Kleinig TJ, *et al*. Endovascular therapy for ischemic stroke with perfusion imaging selection. *N Engl J Med* 2015;372:1009–18.
- 4 Saver JL, Goyal M, Bonafé A, *et al*. Stent retriever thrombectomy after intravenous t-PA vs. t-PA alone in stroke. *N Engl J Med* 2015;372:2285–95.
- 5 Jovin TG, Chamorro A, Cobo E, *et al*. Thrombectomy within 8 hours after symptoms onset in ischemic stroke. *N Engl J Med* 2015;372:2296–306.
- 6 Kriscsik O, Miloslavski E, Fischer S, *et al*. A comparison of functional and physical properties of self-expanding intracranial stents (Neuroform 3, Wingspan, Solitaire, Leo+, Enterprise). *Minim Invas Neurosurg* 2011;54:21–8.
- 7 Ranc H, Servais C, Chauvy PF, *et al*. Effect of surface structure on frictional behaviour of a tongue/palate tribological system. *Tribology Int* 2006;39:1518–26.
- 8 Boeckh-Behrenc T, Schubert M, Forschler A, *et al*. The impact of histological clot composition in embolic stroke. *Clin Neuroradiol* 2014. doi:10.1007/s00062-014-0347-x
- 9 Chueh JY, Khun AL, Puri AS, *et al*. Reduction in distal emboli with proximal flow control during mechanical thrombectomy. *Stroke* 2013;44:1396–401.
- 10 Machi P, Costalat V, Lobotesis K, *et al*. Solitaire FR thrombectomy system: immediate results in 56 consecutive acute ischemic stroke patients. *J Neurointerv Surg* 2012;4:62–6.
- 11 Turk AS, Frei D, Fiorella D, *et al*. ADAPT FAST Study: a direct aspiration first pass technique for acute stroke thrombectomy. *J Neurointerv Surg* 2104;6:260–4.
- 12 Katori N, Tanaka KA, Szlam F, *et al*. The effect of platelet count on clot retraction and tissue-plasminogen activator-induced fibrinolysis on thrombelastography. *Anesth Analg* 2005;100:1781–5.
- 13 Ogata J, Yutani C, Otsubo R, *et al*. Heart and vessel pathology underlying brain infarction in 142 stroke patients. *Ann Neurol* 2008;63:770–81.
- 14 Wenger KJ, Berkefeld J, Wagner M, *et al*. Flat panel detector computed tomography for the interaction between contrast-enhanced thrombi and stent retrievers in stroke therapy: a pilot study. *Clin Neuroradiol* 2014;24:251–4.
- 15 Wenger KJ, Nagl F, Wagner M, *et al*. Improvement of stent retriever design and efficacy of mechanical thrombectomy in a flow model. *Cardiovasc Intervent Radiol* 2013;36:192–7.
- 16 Madjidyar J, Hermes J, Freitag-Wolf S, *et al*. Stent-thrombus interaction and the influence of the aspiration on mechanical thrombectomy: evaluation of different stent retrievers in a circulation model. *Neuroradiology* 2015;57:791–7.
- 17 Mordasini P, Brekenfeld C, Byrne J, *et al*. Experimental evaluation of immediate recanalization effect and recanalization efficacy of a new thrombus retriever for acute stroke treatment in vivo. *AJNR A J Neuroradiol* 2013;34:153–8.
- 18 Gralla J, Schroth G, Remonda L, *et al*. Mechanical thrombectomy for acute ischemic stroke thrombus-device interaction, efficiency, and complications in vivo. *Stroke* 2006;37:3019–24.



Experimental evaluation of stent retrievers' mechanical properties and effectiveness

Paolo Machi, Franck Jourdan, Dominique Ambard, Cedric Reynaud, Kyriakos Lobotesis, Mathieu Sanchez, Alain Bonafé and Vincent Costalat

J NeuroIntervent Surg 2017 9: 257-263 originally published online March 25, 2016

doi: 10.1136/neurintsurg-2015-012213

Updated information and services can be found at:
<http://jniss.bmj.com/content/9/3/257>

These include:

References

This article cites 17 articles, 3 of which you can access for free at:
<http://jniss.bmj.com/content/9/3/257#BIBL>

Open Access

This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/4.0/>

Email alerting service

Receive free email alerts when new articles cite this article. Sign up in the box at the top right corner of the online article.

Topic Collections

Articles on similar topics can be found in the following collections

[Ischemic stroke](#) (414)
[Open access](#) (73)

Notes

To request permissions go to:
<http://group.bmj.com/group/rights-licensing/permissions>

To order reprints go to:
<http://journals.bmj.com/cgi/reprintform>

To subscribe to BMJ go to:
<http://group.bmj.com/subscribe/>